




Eastern  
**Health**

**EASTERN HEALTH LABORATORY**

**TEST DICTIONARY**

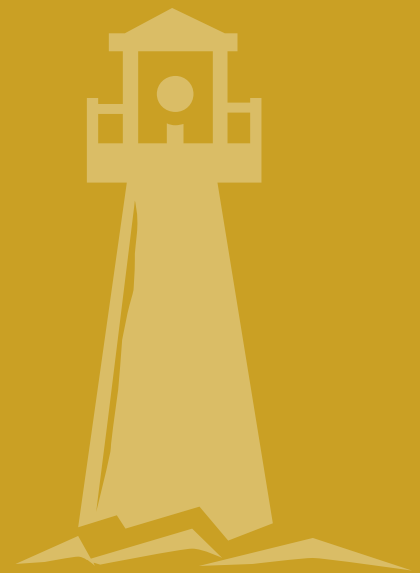

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
1, 25 DIHYDROXY VIT D	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
7 DEHYDROCHOLESTEROL	COLLECT 3.5 ML OF BLOOD IN LITHIUM HEPARIN (LIGHT GREEN) OR COLLECT 1 LITHIUM HEPARIN CONE  PROTECT FROM LIGHT FREEZE PLASMA
11 DESOXYCORTISOL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
17 HYDROXYPROGESTERONE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
25 HYDROXY VIT D	COLLECT 4.0 BLOOD IN SERUM TUBE (RED)  CENTRIFUGE AND REFRIDGERATE SERUM
ACANTHOMOEBA	CORNEALSCRAPING/SWAB IN 1 ML OF STERLE SALINE  SEND CONTACT LENS CASE/FLUID
ACETAMINOPHEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  NOTE: RECORD TIME DOSAGE WAS GIVEN
ACETONE SERUM (QUANT)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ACETYLCHOLINE RECEPTOR AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ACTH STIM.FOR CORT.	COLLECT 3.5 ML OF BLOOD IN LITHIUM TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
ADENOVIRUS	COLLECT RESPIRATORY, CONJUNCTIVAL OR BODY FULIDS IN UTM TRANSPORT MEDIUM
ADRENAL AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
ADRENOCORTICOTROPHIC HORMONE	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) DELIVER TO LAB ON ICE
ALANINE AMINOTRANSFERASE (ALT)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ALBUMIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ALBUMIN (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
ALDOLASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ALDOSTERONE	COLLECT 4.0 ML OF BLOOD IN SST TUBE (GOLD)
ALK PHOS ISOENZYMES	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ALKALINE PHOSPHATASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ALPHA 1 ANTITRYPSIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ALPHA FETOPROTEIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE
ALPHA GLYCOPROTEIN SUBUNIT	COLLECT 3.5 ML OF BLOOD IN SST (GOLD)
ALPHA THALASSEMIA	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) TUBE MUST BE FULL
ALUMINUM	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
AMIKACIN PEAK	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT 30 MINIUTES AFTER END OF INFUSION. RECORD TIME OF COLLECTION

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
AMIKACIN TROUGH	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT JUST PRIOR TO NEXT DOSE
AMINOACIDS PLASMA QUANTITATIVE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN) OR TWO SODIUM HEPARIN CONES PROTECT FROM LIGHT CENTRIFUGE AND FREEZE PLASMA
AMIODARONE (CORDARONE)	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
AMITRIPTYLINE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
AMMONIA	COLLECT 4.0 ML SODIUM HEPARIN TUBE (DARK GREEN) PEDIATRIC CONES ARE ACCEPTABLE. FILL COMPLETEL
AMNIOTIC FLUID BILIRUBIN SCAN	COLLECT 5 ML OF AMNOITIC FLUID IN STERILE CONTAINER PROTECT SPECIMEN FROM LIGHT, KEEP COLD BUT DO NOT FREEZE
AMOXICILLIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
AMPICILLIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
AMYLASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
AMYLASE (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
ANA	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ANCA BY ELISA	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ANDROSTENEDIONE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
ANGIOTENSIN CONVERTING ENZYME	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ANTIBODIES TO TETANUS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ANTI-CARDIAC MUSCLE AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ANTIDIURETIC HORM. (VASOP)	COLLECT 3 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
ANTI-DNA AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ANTI-HBe	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CAN BE PERFORMED ON THE SAME TUBE AS HBSAG
ANTIMULLAERIAN HORMONE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ANTITHROMBIN 111(F)	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
ANTI-THYROID PERIODASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
APO-A1	COLLECT 3.5 ML OF BLOOD IN SST TUBE PATIENT MUST BE FASTING FOR 12 HR
APOLIPOPROTEIN B	. 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROVIDE CLINICAL HISTORY OF PATIENT
ARSENIC	COLLECT 4.0 ML BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
------------------	---------------------

---

ARSENIC URINE TOTAL 24 HOUR	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER  ADJUST 100 ML ALIQUOT PH TO 2-3 WITH 6N HCL AFTER COLLECTION  PATIENT TO AVOID SEAFOOD CONSUMPTION 5 DAYS PRIOR TO TEST
-----------------------------	---

---

ARTERIAL BLOOD GASES	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPILLARY TUBE.  MUST DELIVER TO CHEMISTRY LAB IMMEDIATELY
----------------------	---

---

ARTERIAL CORD GAS	COLLECT A MINIMUM OF 0.5 ML OF BLOOD IN HEPARINIZED SYRINGE
-------------------	---

---

ASCARIS ANTIBODIES	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  PROVIDE CLINICAL HISTORY OF PATIENT
--------------------	---

---

ASPARTATE AMINOTRANSFERASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
----------------------------	--

---

ASPERGILLUS PRECIPITAN AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
---------------------------	--

---

AURAMINE CONCENTRATED SMEAR	COLLECT SPUTUM, URINE, BRONCHIAL WASH/LAVAGE, TISSUE, FLUID IN STERILE CONTAINER  (THIS TEST IS COMPLETED AS PART OF THE AFB/TB CULTURE GROUP TEST.)
-----------------------------	--

---

B-2-GLYCOPROTEIN AB	COLLECT 2 X 3.5 ML OF BLOOD IN SST TUBE (GOLD) TUBES MUST BE FULL
---------------------	---

---

BARTONELLA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
---------------------	--

---

BENZENE	COLLECT 4 X 3.5 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN)  DO NOT FREEZE
---------	---

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
BETA 2-MICROGLOBULIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
BETA HUMAN CHORIONIC GONAD	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
BETA HYDROXY BUTYRATE	COLLECT 3.5 ML OF BLOOD SST TUBE (GOLD)  PLACE ON ICE IMMEDIATELY CENTRIFUGE, SEPARATE AND FREEZE UNTIL PROCESSING
BETA THALASSEMIA	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) TUBES MUST BE FULL
BG+800 PANEL BLOOD	COLLECT ONE FULL CAPILLARY TUBE OF BLOOD. DELIVER TO LAB IMMEDIATELY
BILE ACIDS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  CENTRIFUGE, SEPARATE AND FREEZE SERUM
BILIRUBIN PROFILE NEONATAL	COLLECT 600 UL OF BLOOD IN LIGHT GREEN CONE
BIOAVAILABLE TESTOSTERONE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
BIOTINIDASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
BLASTOMYCES SEROLOG Y	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
BLOOD BILIRUBIN	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPILLARY TUBE (LIGHT GREEN)  MUST DELIVER TO CHEMISTRY LAB IMMEDIATELY
BLOOD ETHANOL	COLLECT 3.5 ML OF BLOOD IN SODIUM FLORIDE TUBE (GREY)  TO BE USED FOR FORENSIC SPECIMENS ONLY

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
BLOOD FILM HAEMATOLOGIST	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) COLLECT BEFORE 1000 HRS MONDAY TO WEDNESDAY
BLOOD FILM MORPHOLOGY	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
BLOOD GLUCOSE 2 HR	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPELLARY TUBE MUST DELIVER TO LAB IMMEDIATELY
BLOOD GLUCOSE/LACTATE	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPELLARY TUBE MUST DELIVER TO LAB IMMEDIATELY
BLOOD LACTATE	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPELLARY TUBE
BODY FLUID BARB.SCREENBF	COLLECT 20 ML OF FLUID IN STERILE URINE CONTAINER
BODY FLUID COUNT & DIFF	PLEASE CALL HAEMATOLOGY LAB FOR INSTRUCTIONS 777-6386
BODY FLUID MAGNESIUM	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN) MUST DELIVER TO LAB IMMEDIATELY
BONE MARROW DIFF	PLEASE CALL HAEMATOLOGY LAB FOR INSTRUCTIONS 777-6386
BONE MARROW DONOR	THIS IS MICROBIOLOGY REQUEST CALL MICROBIOLOGY AT 777-6936
BREATH HYDROGEN TEST	THIS TEST MUST BE BOOKED BY MASS SPECTROMETRY LAB AT 777- 4511 OR 777-4513



**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
BRUCELLA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CALL PUBLIC HEALTH LAB IMMEDIATLEY AT 777-6536
BUFFY COAT SMEARS	COLLECT 1 ML OF BLOOD IN EDTA TUBE (LAVENDER)
C DIFFICILE	COLLECT 5 ML OF FLUID FECES OR 3 G FORMED STOOL IN STERILE CONTAINER NO PRESERVATIVE
C.DIFFICILE TOXIN A/B	COLLECT 5 ML OF FLUID FECES OR 3 G FORMED STOOL IN STERILE CONTAINER
C.DIPHThERIA ANTITOXIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
C1 EST. INH. FUNCTIONAL ASSAY	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) ON ICE
C1 ESTERASE INHIBITOR	COLLECT 4.5 ML OF PLASMA IN NA CITRATE TUBE (BLUE) DELIVER TO LAB ON ICE CENTRIFUGE, SEPARATE AND FREEZE
CA 15-3	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CA19-9	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
CADMIUM	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN) REFRIGERATE BLOOD
CAFFEINE	COLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CALCITONIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) DELIVER TO LAB ON ICE - PATIENT MUST BE FASTING
CALCIUM	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
CALCIUM (IONIZED)	COLLECT 3.5 ML OF BLOOD IN SST TUBE DO NOT USE TOURNIQUET - DELIVER TO LAB ON ICE
CALCIUM (IONIZED) BLOOD	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) DO NOT USE TOURNIQUET AND DELIVER ON ICE TO LAB IMMEDIATELY
CANNABINOID MONITORING	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
CANNABINOIDS BY GCMS	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
CAPILLARY BLOOD GAS	COLLECT ONE FULL CAPILLARY TUBE OF BLOOD DELIVER TO LAB IMMEDIATELY
CARBAMAZEPINE (TEGRETOL)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT PRIOR TO DOSAGE NOTE: RECORD TIME OF COLLECTION
CARBAMAZEPINE EPOXIDE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) DELIVER TO CHEMISTRY LAB A.S.A.P
CARBOHYDRATE DEFIC TRANSFERRIN	COLLECT 4.0 ML OF BLOOD IN (RED TUBE) ONLY
CARBOXYHEMOGLOBIN BY BG	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
CARCINOEMBRYONIC ANTIGEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE
CARDIAC SCREEN (EMERG)	COLLECT 3.5 ML OF BLOOD IN SST (GOLD)
CAROTENE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROTECT FROM LIGHT. AVOID HEMOLYSIS

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
CATECHOLAMINE UR FRACTIONATION	URINE TO BE COLLECTED IN A 24 HR. PRE-ACIDIFIED CONTAINER (50% ACETIC ACID)  NOTE: PLEASE STATE VOLUME AND TIME OF COLLECTION ON CONTAINER
CATECHOLAMINES PLASMA	BY APPOINTMENT ONLY THROUGH CLINIC INVESTIGATION UNIT AT 777-7270
CBC	COLLECT 4.0 ML OF BLOOD IN EDTA (LAVENDER)
CBC / RETIC	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
CD3 COUNT BY FLOW CYTOMETRY	CALL FLOW CYTOMETRY LAB FOR INSTRUCTIONS AT 777-7704
CD34 COUNT BY FLOW CYTOMETRY	CALL FLOW CYTOMETRY LAB FOR INSTRUCTIONS AT 777-7704
CD4/CD8 RATIO	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) CBC MUST ALSO BE PERFORMED AT SAME TIME  TESTING PERFORMED MONDAY-THURSDAY ONLY (EXCEPT BY REQUEST) DELIVER TO LAB ASAP  SAMPLE MUST BE RECEIVED WITHIN 30 HOURS OF COLLECTION.
CERULOPLASMIN	COLLECT 3.5 ML OF BLOOD IN SST (GOLD)
CHLAMYDIA CULTURE	RESPIRATORY OR CONJUNCTIVAL SAMPLE IN UTM TRANSPORT MEDIUM
CHLAMYDIA GROUP SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CHLAMYDIA PSITTACI SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CHLAMYDIA TRACHOMATIS DNA	COLLECT ENDOCERVICAL OR URINE IN COBAS PCR TRANSPORT MEDIUM  INSTRUCTIONS INCLUDED WITH TRANSPORT KIT.

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
CHLORDIAZEPOXIDE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CHLORIDE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CHLORIDE (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
CHLORPROMAZINE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CHOLESTEROL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) 12 HOUR FAST IS REQUIRED
CHOLESTEROL (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
CHROMIUM	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN)
CHROMIUM URINE	COLLECT IN 24 HR URINE CONTAINER STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER. REFRIGERATE
CHROMOGRANIN A	COLLECT 4.0 ML OF BLOOD IN (RED) TUBE DO NOT COLLECT IN SST TUBE
CITALOPRAM (CELEXA)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CLOBAZAM (FRISIUM)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CLOMIPRAMINE (ANAFRANIL)	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
CLOSTRIDIUM BOTULINUM	PLEASE CALL PUBLIC HEALTH LAB FOR INSTRUCTIONS AT 777-6536

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
CLOSTRIDIUM BOTULINUM SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CALL PUBLIC HEALTH LAB IMMEDIATLEY AT 777-6536
CLOZAPINE (CLOZARIL)	COLLECT 3.0 ML OF BLOOD IN LITHIUM HEPARIUM TUBE (LIGHT GREEN)
CMV CULTURE	COLLECT URINE, RESPIRATORY SECRETIONS, TISSUES, ASPIRIATES, AND OTHER BODY FLUIDS IN UTM TRANSPORT MEDIUM.
CMV IgG	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CMV IgM	COLLECT 4.0 ML OF BLOOD IN RED TUBE
CMV PCR	COLLECT CSF IN STERILE CONTAINER DO NOT ORDER ON BLOOD ORDER CMVL (CMV VIRAL LOAD)
CMV VIRAL LOAD	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) DO NOT ORDER CMV PCR
CO2	COLLECT 3.5 ML OF IN SST TUBE (GOLD) SAMPLE SHOULD BE COLLECTED ANAEROBICALLY
CO2 (FLUID)	COLLECT 3.5 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
COARSE GRAN.CAST	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
COBALT	COLLECT 3.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
COBALT RANDOM URINE	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER, REFRIGERATE

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
COBALT URINE 24 HR.	COLLECT IN 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME ON CONTAINER REFRIGERATE & DELIVER TO LAB
COLD AGGLUTININS	COLLECT 4.0 ML OF BLOOD IN (RED) TUBE  KEEP SAMPLE WARM (37°C) BRING TO LAB IMMEDIATELY
COMB.PIT.STIM.TEST	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)  DELIVER TO LAB ON ICE
COMPLEMENT C2	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  DELIVER TO LAB A.S.A.P
COMPLEMENT C3	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
COMPLEMENT C4	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CONF TYPE IIB VWF	ARRANGE WITH COGULATION LAB FOR TESTING INSTRUCTIONS AT 777-7178
COPPER	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
CORTISOL	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
CORTISOL 0800 HR	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)  THIS SAMPLE MUST BE COLLECTED 0700-0900 OR CANNOT BE TESTED
CORTISOL PM	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
C-PEPTIDE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) DELIVER TO LAB ON ICE ASAP. PATIENT MUST BE FASTING
C-REACTIVE PROTEIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CREATINE PHOSPHOKINASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CREATININE (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
CREATININE CLEARANCE	COLLECT 24 HR. URINE AND 1 3.5 ML SST TUBE (GOLD) STATE VOLUME, DATE AND TIME OF COLLECTION ON URINE CONTAINER
CREATININE/EST. GFR BY MDRD	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CRP HIGH SENSITIVITY	COLLECT 3.5 ML BLOOD IN SST TUBE (GOLD)
CRYOGLOBULIN	COLLECT 4.0 ML OF BLOOD IN RED STOPPER TUBE KEEP SPECIMEN AT 37°C AND DELIVER TO LAB IMMEDIATELY
CRYSTAL, FLUID	COLLECT 4.0 ML OF FLUID IN SODIUM HEPARIN TUBE (DARK GREEN) OR IN PLAIN RED STOPPER TUBE
CSF CHLORIDE	COLLECT 0.5 M L FLUID IN # 2 CSF TUBE
CSF DIFFERENTIAL	COLLECT 0.5 M L OF FLUID IN # 4 CSF TUBE
CSF ERYTHROCYTE COUNT	COLLECT 0.5 M L OF FLUID IN # 4 CSF TUBE
CSF GLUCOSE	COLLECT 0.5 MLOF FLUID IN # 2 CSF TUBE
CSF LACTIC ACID	COLLECT 0.5 M LOF FLUID IN # 2 CSF TUBE

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
CSF LEUKOCYTE COUNT	COLLECT 0.5 M L OF FLUID IN # 4 CSF TUBE
CSF OLIGOCLONAL BANDING	COLLECT 0.5 M LOF FLUID IN # 2 CSF TUBE
CSF PROTEIN	COLLECT 0.5 M LOF FLUID IN # 2 CSF TUBE
CYANIDE	COLLECT 4.0 ML OF BLOOD IN EDTA (LAVENDER)
CYCLIC CITRULLINATED PEPTID AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
CYCLOSPORIN A	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
CYSTINE IN WBC	NOTIFY MASS SPECTROMETRY LAB ONE WEEK IN ADVANCE AT 777-7513 WHEN ORDERING TEST
CYTOTOXIC ABS – HALIFAX	COLLECT 4.0 ML OF BLOOD IN SERUM TUBE (RED)
D2 (DUST MITE)	COLLECT 3.5 ML BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
D-DIMER	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
DEOXYPYRIDINOLINE	COLLECT 20 ML OF FIRST OR SECOND MORNING URINE PUT IN STERILE URINE CONTAINER KEEP COLD
DESIPRAMINE	COLLECT 3.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
DG TYPING BY DNA	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
DHEA SULPHATE	COLLECT 3.5 ML OF BLOOD SST TUBE (GOLD)
DIAZEPAM	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)



**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
DIBUCAINE NUMBER	COLLECT 3.5 ML OF BLOOD IN SST (GOLD)
DIC	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
DIGOXIN	COLLECT 4.0 ML OF BLOOD IN RED STOPPER TUBE DO NOT USE SST TUBE    RECORD TIME OF COLLECTION
DIHYDROPTERIDINE REDUCTASE	NOTIFY MASS SPECTROMETRY LAB IN ADVANCE AT 777-4513 WHEN ORDERING TEST
DIHYDROTESTOSTERONE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
DIRECT BILIRUBIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
DISOPYRAMIDE GREEN)	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
DONATH-LANDSTEINER	COLLECT 4.0 ML OF BLOOD IN (RED) TUBE
DOXEPIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
DRUG SCREEN (URINE) CONTAINER	COLLECT 20 ML OF RANDOM URINE IN STERILE URINE CONTAINER
E. COLI 0157	COLLECT STOOL IN STERILE DRY CONTAINER
E. COLI 0157 VEROTOXIN	COLLECT STOOL IN STERILE DRY CONTAINER
E1 (CAT)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
E2 (DOG)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
EBV PCR	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) COLLECT CSF IN STERILE CONTAINER. DELIVER TO PUBLIC HEALTH LAB WITHIN 24 HOURS
EBV VIRAL LOAD	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
ECHINOCOCCUS SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ELECTROLYTES	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
EMERG.TOXIC SCREEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ENTAMOEBA HYSTOLYTICA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROVIDE CLINICAL HISTORY OF PATIENT
ENTEROVIRUS	CSF OR STOOL IN STERILE CONTAINER THROAT SWAB, BODY FLUID OR TISSUE IN UTM TRANSPORT MEDIUM
EOSINOPHILS NASAL SMEAR	COLLECT SAMPLE IN FRESH NASAL SWAB
EOSINOPHILS SPUTUM	SEND FRESH SPUTUM SPECIMEN IN STERILE CONTAINER
EOSINOPHILS STOOL SMEAR ONLY	COLLECT 1.0 G IN STERILE CONTAINER SEND FRESH SWAB
EOSINOPHILS URINE CONTAINER	COLLECT 20 ML OF FRESH URINE IN STERILE URINE
EPSTEIN-BARR AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
ERY.PROTOPORPHYRIN	COLLECT 4.0 ML BLOOD IN EDTA TUBE (LAVENDER) PROTECT FROM LIGHT. ORDER HEMATOCRIT
ESTRADIOL	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
ESTRONE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ETHANOL PLASMA OR SERUM	COLLECT 3.5 ML OF BLOOD SST TUBE (GOLD) NONALCOHOLIC GERMICAL SOLUTION SHOULD BE USED TO SWAB VENIPUNCTURE SITE
ETHOSUXIMIDE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ETHYLENE GLYCOL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
EUGLOBULIN LYSIS TIME	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) PLACE ON ICE AND DELIVER TO COAGULATION LAB ASAP
E XTRACTABLE NUCLEAR AG.	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
F1 (EGG WHITE)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
F13 (PEANUT)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
F14 (SOYA)	COLLECT 3.5 ML OF WHOLE BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
F17 (HAZELNUT)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

**EASTERN HEALTH TEST DICTIONARY****TEST NAME****INSTRUCTIONS**

F2 (MILK )

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

F20 (ALMOND)

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

F23 (CRAB)

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

F24 (SHRIMP)

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

F3 (FISH)

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

F4 (WHEAT)

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

F5 (RYE)

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

F6 (BARLEY)

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

F7 (OAT)

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

F80 (LOBSTER)

COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
F9 (RICE)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
FACTOR II ASSAY	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR IX ASSAY	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR IX INHIBITOR	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR V ASSAY	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR V INHIBITOR	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR VII ASSAY	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL.
FACTOR VII INHIBITOR	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR VIII COAGULANT	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR VIII INHIBITOR	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR VIII MULTIMERS	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR X ASSAY	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR XI ASSAY	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
FACTOR XI INHIBITOR	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR XII ASSAY	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR XIII ANTIGEN	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR XIII CS	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FACTOR XIII INHIBITOR	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FASTING BLOOD GLUCOSE	COLLECT 3.5 ML OF BLOOD IN SST TUBE OR 4.0 ML SODIUM FLUORIDE POTASSIUM OXALATE TUBE (GREY)  PATIENT MUST BE FASTING.
FECAL ELASTASE	COLLECT 5 GRAMS OF FORMED STOOL IN STERILE CONTAINER
FECAL FAT 72 HOUR	OBTAIN PRE-WEIGHED CONTAINER FROM LAB  72 HOUR COLLECTION. CALL 777-6398
FECAL IMMUNOASSAY TEST	SPECIAL FECAL TUBE (CONTANING NA AZIDE) PROVIDED BY LAB
FECAL PORPHYRINS	COLLECT RANDOM STOOL IN STERILE CONTAINER  PROTECT FROM LIGHT
FETAL FIBRONECTIN	SPECMIN TRANSPORT TUBE MUST CONTAIN SWAB TIP  MUST BE LABELED.

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
FETAL HEMOGLOBIN BY BG	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPILLARY TUBE  DELIVER TO LAB IMMEDIATELY
FETAL HGB SCREEN TEST	COLLECT 4.0 ML OF STOOL IN EDTA TUBE (LAVENDER)
FETAL LUNG MATURITY	COLLECT 10 ML OF AMNIOTIC FLUID IN STERILE CONTAINER.
FIBRINOGEN-CLAUSS	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FILARIA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  PROVIDE CLINICAL HISTORY OF PATIENT
FK-506 (TACROLIMUS)	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) OR ONE LAVENDER CONE
FLETCHER FACTOR	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
FLUID CALCIUM	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR RED STOPPER TUBE
FLUID ETHANOL	COLLECT 3.0 ML FLUID OF FLUID IN LITHIUM HEPARIN (LIGHT GREEN)
FLUORIDE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
FLUOXETINE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN)  COLLECT PRIOR TO AM DOSE
FOLLICLE STIMUL.HORMONE	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
FREE CARBAMAZEPINE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
------------------	---------------------

---

FREE PHENYTOIN	COLLECT 4.0 ML OF BLOOD IN RED TUBE. DO NOT USE SST TUBE  RECORD TIME OF COLLECTION
----------------	---

---

FREE PSA	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
----------	--

---

FREE T3	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
---------	--

---

FREE THYROXINE	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
----------------	--

---

FRUCTOSAMINE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
--------------	--

---

FUNGAL CULTURE	COLLECT HAIR, SKIN OR NAILS IN FUNGAL COLLECTION KIT PROVIDED BY PUBLIC HEALTH LAB
----------------	--

---

GABAPENTIN (NEURONTIN)	COLLECT 4.0 ML OF BLOOD IN SERUM TUBE (RED)  DO NOT USE GEL-SEPARATOR TUBE
------------------------	--

---

GALACTOCEREBROSIDASE	COLLECT TWO SODIUM HEPARIN TUBES (DARK GREEN) ONLY COLLECT FROM 8:00 AM TO 1:30  CALL 777-4513 FOR INSTRUCTIONS
----------------------	--

---

GALACTOSEMIA SCREEN	COLLECT ON FILTER PAPER; LET DRY THEN FREEZE AFTER HOURS  DO NOT COLLECT ON FRIDAYS OR WEEKENDS
---------------------	---

---

GAMMA GLUT TRANSFA	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
--------------------	--

---

GASTRIC CONTENT TOXICOLOGY	COLLECT GASTIC CONTENT IN STERILE CONTANER
----------------------------	--

---

GASTRIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  DELIVER TO LAB ON ICE. PATIENT MUST BE FASTING
---------	--



**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
GENTAMICIN TROUGH	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  COLLECT PRIOR TO NEXT DOSE NOTE: RECORD TIME OF COLLECTION
GESTATIONAL GLUCOSE TOLERANCE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) OR SODIUM FLORIDE TUBE (GREY)  PLEASE USE 100 gm. DRINK FOR THIS TEST  PLEASE RECORD TIMES OF COLLECTION ON TUBE
GLOMERULAR BASEMENT. MEMBRANE AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
GLUCAGON	PRECHILL TUBES, COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)  DELIVER TO LAB ON ICE
GLUCOSE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
GLUCOSE 6 PHOS DEHYDROGENASE	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
GLUCOSE FLUID	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
GLUCOSE TOLERANCE TEST	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  PATIENT TAKES 75 gm. DRINK  PLEASE RECORD TIMES OF COLLECTION ON TUBE
GLUTAMIC ACID DECARBOXYLASE AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
GLYBURIDE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
GRANTS TEST	COLLECT 2 X 8.5 ML OF BLOOD IN ACD TUBE (YELLOW)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
GROWTH HORM.REL.HORMONE.TEST	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) DELIVER TO LAB ON ICE
GROWTH HORMONE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) DELIVER TO LAB ON ICE
H2 (HOUSE DUST)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HAEM COMP-CLASSICAL PATH	COLLECT 3.5 ML OF BLOOD IN (RED) TUBE MUST BE DELIVERED TO LAB ON ICE WITHIN 1 HOUR
HAEMOGLOBIN F	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) TUBE MUST BE FULL
HAEMOGLOBULIN CAST	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
HANTAVIRUS SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROVIDE CLINICAL HISTORY OF PATEINT. ALERT PUBLIC HEALTH LAB
HAPTOGLOBIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HBV VIRAL LOAD	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
HCV VIRAL LOAD	COLLECT 2 X 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HDL CHOLESTEROL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) 12 HR. FAST OR NON-FASTING ACCEPTABLE IF NOT PART OF LIPP PROFILE
HEAT SHOCK PROTEIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HEINZ BODIES	COLLECT 4 ML OF BLOOD IN EDTA TUBE (LAVENDER)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
HELICOBACTER PYLORI IgG	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HEMAGLOBIN A1C	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
HEMOGLOBIN BY BG	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPILLARY TUBE.  DELIVER TO LAB IMMEDIATELY
HEMOGLOBIN BY BLOOD GAS METHOD	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPILLARY TUBE  DELIVER TO LAB IMMEDIATELY
HEMOGLOBIN ELECTRO GROUP	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA (LAVENDER TUBE MUST BE FULL)
HEMOGLOBIN S	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA (LAVENDER) TUBE MUST BE FULL
HEMOGLOBINOPATHY CONFIRMATION	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) TUBE MUST BE FULL
HEMOPHILIA STUDY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HEMOPHILUS INFLUENZA B AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HEMOSIDERIN SPUTUM	COLLECT 0.5 ML OF SPUTUM IN STERILE SWAB OR STERILE CONTAINER
HEMOSIDERIN URINE	COLLECT 20 ML OF FRESH URINE IN STERILE URINE CONTAINER
HEP IND THROMBO COMFIRMATION	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) TUBE MUST BE FULL
HEP IND THROMBOCYTOPENIA	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) TUBE MUST BE FULL

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
HEPATITIS A AB (IgG)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CAN BE PERFORMED ON SAME TUBE AS HBSAG
HEPATITIS A AB (IgM)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CAN BE PERFORMED ON THE SAME TUBE AS HBSAG
HEPATITIS B CORE AB-TOTAL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CAN BE PERFORMED ON SAME SAMPLE AS HBSAG
HEPATITIS B CORE AB-IgM	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CAN BE PERFORMED ON SAME TUBE AS HBSAG
HEPATITIS B DNA PCR	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
HEPATITIS B E" Ag"	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CAN BE PERFORMED ON THE SAME TUBE AS HBSAG
HEPATITIS B GENOTYPING	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
HEPATITIS B SURFACE AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CAN BE PERFORMED ON SAME TUBE AS HBSAG
HEPATITIS B SURFACE ANTIGEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HEPATITIS C "RNA (PCR)"	COLLECT 3 X 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT 2 X 3.5 ML OF BLOOD IN SST TUBE (GOLD) CAN BE PERFORMED ON SAME TUBE AS HBSAG
HEPATITIS C GENOTYPING	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
HEPATITIS C LIA	COLLECT 3.5 ML OF BLOOD IN SST TUBE CAN BE PERFORMED ON SAME TUBE AS HBSAG
HEPATITIS D (ANTI-DELTA)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) CAN BE PERFORMED ON THE SAME TUBE AS HBSAG
HEPATITIS E SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HEPATITIS SCREEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HERPES SIMPLEX ANTIGEN	COLLECT VESICULAR FLUID, SCRAPINGS FROM BASE OF SKIN LESIONS USING MICROTRACK HSV KIT ACCOMPANIED BY SAMPLE FOR HSVCR (CULTURE)
HERPES SIMPLEX CULTURE	COLLECT VESICULAR FLUID, SCRAPINGS FROM BASE OF SKIN LESIONS, BIPOSY MATERIAL IN UTM TRANSPORT MEDIUM.
HEXOSAMINIDASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HIGH MOLECULAR WEIGHT KINNOGEN	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE)
HISTAMINE PLASMA	PRE-CHILLED TUBE, COLLECT 4 ML OF BLOOD IN EDTA TUBE (LAVENDER) DELIVER TO LAB ON ICE
HISTAMINE URINE	COLLECT URINE IN A 24 HR URINE CONTAINER STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
HISTONE AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HISTOPLASMA AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
HIV 1 AB IMMUNOBLOT	PERFORMED ON SAME SAMPLE AS HIV ANTIBODY
HIV AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HIV GENOTYPE	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) SPECIMEN MUST BE COLLECTED AS EARLY AS POSSIBLE IN DAY BEFORE 12:30 PM * SEND TO MICRO HSC STAT
HIV PCR	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) CALL PUBLIC HEALTH LAB BEFORE TAKING TEST. SAMPLE MUST BE SENT TO OTTAWA SAME DAY.
HIV VIRAL LOAD	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) SPECIMEN MUST BE COLLECT AS EARLY AS POSSIBLE IN THE DAY, BEFORE 12:30 PM * SEND TO MICRO. HSC STAT
HLA B27	COLLECT 2 X 4.0ML OF BLOOD IN EDTA TUBE (LAVENDER)
HOMOCYSTEINE	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) DELIVER TO LAB SORTING.CENTRIFUGE ASAP AND FREEZE PLASMA
HOMOVANILLIC ACID	COLLECT URINE IN A 24 HR URINE CONTAINER STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
HSV PCR	COLLECT CSF IN STERILE CONTAINER
HSV SEROLOGY DONOR SCREEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HSV SEROLOGY GENITAL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
HTLV I & II	COLLECT 3.5 OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
HTLV I & II PCR	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)  CALL PUBLIC HEALTH LAB BEFORE TAKING TEST. SAMPLE MUST BE SENT TO OTTAWA SAME DAY
HYDROXYPROLINE TOTAL	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
HYPERCOAGULABLE SCREEN	COLLECT 3 X 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) AND 1 X 3.5 ML OF BLOOD IN SST TUBE (GOLD). TUBES MUST BE FULL
IBUPROFEN (MOTRIN)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ICTOTEST	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
IGE ALLERGY PANEL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
IGFBP3	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  DELIVER TO LAB ON ICE
IGG SUBCLASSES	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
IMIPRAMINE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
IMMUNOGLOBULIN A	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
IMMUNOGLOBULIN E	COLLECT 3.5 ML OF BLOOD IN SST TUBE(GOLD)
IMMUNOGLOBULIN G	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
IMMUNOGLOBULIN M	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
IMMUNOGLOBULINS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
IN.LIKE GROWTH FACTOR 1	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) DELIVER TO LAB ON ICE
INFLUENZA VIRUS	NASOPHARYNGEAL SWAB OR ASPIRATE IN UTM TRANSPORT MEDIUM TRACHEAL SECRETIONS, BRONCHIAL WASH/LAVAGE OR TISSUE IN STERILE CONTAINER
INHALANTS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
INORGANIC PHOSPHOROUS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) FASTING SAMPLE IS PREFERRED
INSULIN	COLLECT 3.5 ML BLOOD IN SST TUBE (GOLD)
INSULIN AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
INSULIN TOLERANCE TEST	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PLEASE RECORD TIMES OF COLLECTION ON TUBE
INTRINSIC FACTOR AB	COLLECT 3.5 ML BLOOD IN SST TUBE (GOLD)
IRON & IBC CAPACITY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) FASTING A.M. SPECIMEN COLLECTED BEFORE PATIENT IS GIVEN THERAPEUTIC IRON. IRON ANALYSIS SHOULD BE DELAYED SEVEN DAYS IF A TRANSFUSION HAS BEEN GIVEN
ISLET CELL AB.	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ISOPROPANOL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)



**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
------------------	---------------------

---

J.METHYLENEDIOXYMETHAMPH. COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER

---

K82 (LATEX) COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  
NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

---

KLEIHAUER BELKE FETAL HGB COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)

---

LACTATE COLLECT 4.0 ML BLOOD IN SODIUM FLUORIDE TUBE (GREY)  
PLACE ON ICE & DELIVER TO LAB IMMEDIATELY

---

LACTATE DEHYDROGENASE COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  
NO TOURNIQUET. DELIVER TO LAB ASAP

---

LACTIC ACID COLLECT 4.0 ML OF BLOOD IN SODIUM FLUORIDE TUBE (GREY TUBE)  
NO TOURNIQUET. DELIVER TO LAB ON ICE ASAP

---

LAMOTRIGINE (LAMICTAL) COLLECT 4.0 ML OF BLOOD IN SERUM TUBE (RED)

---

LBM AB COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

---

LD (FLUID) COLLECT 4.0 ML OF FLUID IN SODIUM HEPARIN TUBE (DARK GREEN) OR 4.0 ML IN RED STOPPER TUBE

---

LD ISOENZYMES COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

---

LEAD COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)

---

LEGIONELLA ANTIGEN COLLECT 20 ML URINE IN STERILE CONTAINER

---

LEISHMANIA SEROLOGY COLLECT 3.5 ML OF BLOOD IN SST TUBE. (GOLD)  
PROVIDE CLINICAL HISTORY OF PATIENT

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
------------------	---------------------

---

LEPTIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE. (GOLD)  LET CLOT ON ICE FOR 30 MIN. FASTING SPECIMEN. NOTIFY LAB 777-7750
--------	---

---

LEPTOSPIRA	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
------------	--

---

LEUCAEMIA STUDY	COLLECT 3 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)  CALL FLOW CYTOMETRY LAB FOR COLLECTION INSTRUCTIONS AT 777-7704. DELIVER ASAP
-----------------	---

---

LEUK. ALPHA GALACTOSIDASE	COLLECT TWO SODIUM HEPARIN TUBES (DARK GREEN)  ONLY COLLECT 8:00 AM TO 1:30 PM CALL LAB AT 777-4513 DELIVER ASAP
---------------------------	---

---

LEUKOCYTE ARYL SULFATASE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN) ONLY COLLECT BETWEEN 8:00AM -1:30 PM  MONDAY TO FRIDAY. DELIVER TO LAB ASAP. CALL 777-4513
--------------------------	---

---

LEUKOCYTE GALACTOSIDASE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)  ONLY COLLECT BETWEEN 8:00 AM -1:30 PM MONDAY TO FRIDAY DELIVER TO LAB ASAP. CALL 777- 4513
-------------------------	---

---

LEUKOCYTE HEXOSAMINIDASE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)  ONLY COLLECT BETWEEN 8:00 AM - 1:30 PM MONDAY TO FRIDAY DELIVER TO LAB ASAP CALL 777- 4513
--------------------------	---

---

LEUTIC SEROLOGY CSF	COLLECT IN #2 CSF TUBE. PATIENT MUST BE SERUM SYPHILIS POSITIVE
---------------------	--

---

LEVETIRACETAM (KEPPRA)	COLLECT 4.0 ML OF BLOOD IN SERUM TUBE (RED)
------------------------	---

---

LIP.PROF-CHOL, TG, LDL, HDL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  LIPID PROFILIE REQUIRES 12 HR FAST
-----------------------------	--

---

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
LIPASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
LIPOPROTEIN (a)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
LIPOPROTEIN ELECTROPHORESIS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
LITHIUM	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT AT LEAST 12 HR AFTER LAST DOSE
LIVER FUNCTION TESTS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
LONG CHAIN FATTY ACIDS	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
LUPUS ANTICOAGULANT	COLLECT 4.5 ML OF BLOOD IN NA CITRATE (BLUE) TUBE MUST BE FULL
LUTEINIZING HORMONE.	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
LYME IgM WESTERN BLOT	PERFORMED ON SAME SAMPLE AS LYME SEROLOGY IgG/IgM
LYME PCR	CALL PUBLIC HEALTH LAB BEFORE TAKING TEST AT 777-6582
LYME SEROLOGY IgG/IgM	COLLECT 3.5 ML OF BLOOD INN SST TUBE (GOLD)
LYMPHOCYTE CULTURE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)  TESTING PERFORMED MONDAY-THURSDAY ONLY. DELIVER TO LAB ASAP
LYMPHOCYTIC CHORIOMENINGITIS V	SUBMIT SERUM FOR SEROLOGY: CSF OR WHOLE BLOOD IN EDTA FOR PCR  NOTE: RELEVANT CLINICAL HISTORY MUST BE PROVIDED

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
LYMPHOMA STUDY	CALL LAB FOR COLLECTION INSTRUCTIONS AT 777-7704. DELIVER ASAP
MAG AUTOANTIBODY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
MAGNESIUM	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
MALARIA	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
MALARIA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
MANGANESE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
MANGANESE BLOOD	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
MANUAL DIFF	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) REASON FOR ORDERING MUST BE GIVEN
MATERNAL SERUM SCREEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  MUST BE COLLECTED ON PREGNANT PATIENTS THAT ARE BETWEEN 15 WEEKS TO 20 WEEKS 6 DAY GESTATION
MEASLES AB IgM	PERFORMED ON THE SAME SAMPLE AS MEASLES IMMUNE STATUS
MEASLES IMMUNE STATUS (IgG)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
MERCURY	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM (DARK GREEN)
MERCURY RANDOM URINE	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
MERCURY URINE 24 HOUR	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER  ADJUST 100 ML ALIQUOT PH TO 2-3 WITH 6N HCL AFTER COLLECTION BEFORE REFERRAL
METANEPHRINES PLASMA	COLLECT 3 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
METANEPHRINES URINE	URINE TO BE COLLECTED IN A 24 HR. PRE-ACIDIFIED CONTAINER (50% ACETIC ACID)  NOTE: PLEASE STATE VOLUME AND TIME OF COLLECTION ON CONTAINER
METAPNEUMORVIS	NASOPHARYNGEAL SWAB OR ASPIRATE IN UTM TRANSPORT MEDIUM
METHANOL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
METHEMOGLOBIN BY BG	COLLECT 4.0 ML BLOOD IN SODIUM FLUORIDE TUBE (GREY)  PLACE ON ICE & DELIVER TO BIOCHEMISTRY LAB IMMEDIATELY
METHOTREXATE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  PROTECT FROM LIGHT
METHYLMALONIC ACID PLASMA	COLLECT 4.0 ML BLOOD IN EDTA TUBE (LAVENDER) CENTRIFUGE FREEZE PLASMA
MICROSCOPIC	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
MITOCHONDRIAL AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
MIXED GRASS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)  NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
MIXED MOLD	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
MRSA	CALL MICROBIOLOGY LAB FOR INSTRUCTIONS AT 777-6936
MUCOPOLYSACCHARIDES	COLLECT 5 ML OF RANDOM URINE IN NON-STERILE URINE CONTAINER
MUMPS (IgM)	PERFORMED ON THE SAME SAMPLE AS MEASLES IMMUNE STATUS
MUMPS CULTURE	MUMPCR IS THE RECOMMENED METHOD
MUMPS IMMUNE STATUS	PLEASE COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
MUMPS RNA PCR	COLLECT BUCCAL SWAB IN UTM TRANSPORT MEDIUM
MYCOBACTERIAL BL. CULTURE	COLLECT 10 ML OF BLOOD IN DUPONT ISOLATOR
MYCOPHENOLATE (CELLCEPT)	COLLECT 4.0 ML OF PLASMA IN EDTA TUBE (LAVENDER)
MYCOPLASMA IgM AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
NA & K FOR JANEWAY USE ONLY	COLLECT 600 UL OF BLOOD IN LIGHT GREEN CONE
NA & K ON BLOOD	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPILLARY TUBE (LIGHT GREEN) DELIVER TO LAB IMMEDIATELY
NA, K, ICA, CL (BLOOD)	COLLECT 0.2 ML OF BLOOD IN HEPARINIZED SYRINGE OR COLLECT A FULL CAPILLARY TUBE (LIGHT GREEN) DELIVER TO LAB IMMEDIATELY
N-ACETY PROCAINAMIDE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
NEONATAL ALLOIMMUNE	ARRANGE WITH COAGULATION LAB FOR TESTING INSTRUCTIONS AT 777-7178 BOTH MOM AND DAD NEEDED
NEONATAL SCREENING	COLLECT BLOOD ON FILTER PAPER. ALL FOUR CIRCLES SHOULD BE FILLED.  BABIES MUST BE A MINIMUM OF 24 HOURS OLD. RECOMMENDED AGE IS 48 -72 HOUR OLD. ONCE SAMPLE IS  COLLECTED PROTECT FROM HEAT OR SUNLIGHT FOR A LEAST 1 HOUR BEFORE PLACING IN PLASTIC ENVELOPE
NEONATAL T4	THIS TEST IS REFERRED FROM NEONATAL UNITS ONLY.
NEUTROPHIL AB	COLLECT 8.5 ML OF BLOOD IN ACD TUBE (YELLOW) & COLLECT 4.0 ML OF BLOOD IN (RED) TUBE. TESTING  PERFORMED MONDAY TO WEDNESDAY ONLY MUST BE PRE- BOOKED WITH IMMUNOLOGY LAB
NICKEL, PLASMA	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
NMO-IgG	COLLECT 4.0 ML OF BLOOD IN SERUM TUBE (RED)
NORCLOBAZAM	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
NOROVIRUS	COLLECT STOOL IN STERILE DRY CONTAINER
NORTRIPTYLINE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN)
OCCULT BLOOD FECES	SUBMIT ON SPECIAL HEMOCCULT SLIDES ONLY
OCCUPATIONAL EXPOSURE SOURCE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
OCCUPATIONAL EXPOSURE WORKER	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
OCCUPATIONAL HEALTH SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
OSMOLALITY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
OSMOLALITY (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIM TUBE (LIGHT GREEN)
OSMOTIC FRAGILITY	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN) OR LITHIUM HEPARIUM (LIGHT GREEN)  DO NOT COLLECT ON FRIDAYS
OSTEOCALCIN	COLLECT 3.5 ML BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)  DELIVER TO LAB ON ICE
OXIDATIVE BURST	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDAR )  DELIVER TO FLOW CYTOMETRY LAB ASAP
PANCREATIC POLYPEPTIDE	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)  DELIVER TO LAB ON ICE
PARAINFLUENZA VIRUS	NASOPHARYNGEAL SWAB OR ASPIRATE IN UTM TRANSPORT MEDIUM  TRACHEAL SECRETIONS, BRONCHIAL WASH/LAVAGE OR TISSUE IN STERILE CONTANIER
PARANEOPLASTIC AUTOANTIBODIES	COLLECT 4.0 ML OF BLOOD IN (RED) TUBE
PARATHYROID HORMONE	COLLECT 3.5 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
PARIETAL CELL AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PARVOVIRUS B19 SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)



**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
PARVOVIRUS IgG SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PARVOVIRUS IgM SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PATHOLOGICAL CASTS	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
PED SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PEMPHIGUS SCREEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PENICILLIN G	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
PENICILLIN PROFILE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
PENICILLIN V	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
PENTAGASTRIN STIM.TEST	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) DELIVER ON ICE IMMEDIATELY TO RENAL LAB
PH (FLUID)	COLLECT 4.0 ML OF FLUID IN RED STOPPER TUBE COLLECT ANAEROBICALLY AND DELIVER ON ICE TO LAB IMMEDIATELY
PHENOBARBITAL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT AT TROUGH. NOTE TIME, DOSAGE AND TIME COLLECTED
PHENYLALANINE FILTER PAPER	COLLECT ON NEWBORN SCREENING FILTER PAPER. LET AIR DRY 1 HOUR BEFORE PLACING IN ENVELOPE

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
PHENYTOIN (DILANTIN)	COLLECT 4.0 ML OF BLOOD IN SMALL RED TUBE DO NOT USE SST TUBE  RECORD TIME OF COLLECTION
PHYTANIC ACID	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN)  DELIVER TO LAB ON ICE
PI TYPING	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PLASMA ACYL CARNITINE	COLLECT 3.5 ML OF BLOOD IN SODIUM HEPARIN (DARK GREEN) OR THREE LITHIUM HEPARIN CONES  CENTRIFUGE AND FREEZE PLASMA
PLASMINOGEN	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE). TUBE MUST BE FULL
PLATELET AGGREGATION	COLLECT 4 X 2.7 ML PLASTIC NA CITRATE (BLUE)  BOOK WITH COAGULATION LAB.
PLATELET ALLOIMMUNE	ARRANGE WITH COAGULATION LAB FOR TESTING. 777-7178
PLATELET FUNCTION ANALYSIS	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE)  DO NOT SEND IN PNEUMATIC TUBE** MUST BE TESTED WITHIN FOUR HOURS OF COLLECTION
PLATELET INVESTIGATION	ARRANGE WITH COAGULATION LAB FOR TESTING INSTRUCTIONS AT 777-7178
PNEUMOCOCCUS AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PNH STUDY BY FLOW CYTOMETRY	COLLECT 2 X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)  TESTING PERFORMED MONDAY-THURSDAY ONLY. DELIVER TO IMMUNOLOGY LAB ASAP

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
POLYOMAVIRUSES	COLLECT 4.0 ML BLOOD IN EDTA TUBE (LAVENDER) DELIVER TO PUBLIC HEALTH LAB WITHIN 24 HOURS.
POST TRANSFUSION PURPURA	ARRANGE WITH COAGULATION LAB FOR TESTING INSTRUCTIONS AT 777-7178
POTASSIUM (FLUID)	COLLECT 4.0 ML OF FLUID IN SODIUM HEPARIN TUBE (DARK GREEN) OR 4.0 ML IN RED STOPPER TUBE
PREALBUMIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PREGNANCY TEST	COLLECT 20 ML OF "EARLY" MORNING URINE IN A STERILE URINE CONTAINER
PREGNANEDIOL	COLLECT URINE IN A 24 HR URINE CONTAINER WITH BORIC ACID STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
PROCAINAMIDE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
PROGESTERONE	COLLECT 4.0 ML OF BLOOD IN RED STOPPER TUBE DO NOT USE SST TUBE (GOLD)
PROINSULIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PROLACTIN	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD) DELIVER TO CHEMISTRY LAB ASAP
PROSTATE SPECIFIC ANTIGEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
PROTEIN (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
PROTEIN C (F)	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
PROTEIN S(FREE)	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE TUBE MUST BE FULL
PROTEIN TO CREATININE RATIO UR	COLLECT URINE IN A 24 HR URINE CONTAINER STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
PROTRIPTYLINE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN)
PTERIN ANALYSIS	COLLECT 20 ML OF RANDOM URINE IN STERILE URINE CONTAINER THE URINE MUST BE PLACE IN FOIL COVER TUBE CONTAINING 100 MG ASORBIC ACID IMMEDIATELY AFTER COLLECTION BLOOD MUST BE COLLECTED FILTER PAPER
PTH RELATED PEPTIDE	COLLECT 3.5 ML OF BLOOD IN SERUM TUBE (RED) DO NOT USE SST TUBE
PTI PTT	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
PT-INR	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL MUST BE TESTED WITHIN 24 HOURS OF COLLECTION
PTT	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE). TUBE MUST BE FULL SAMPLE MUST BE SPUN WITHIN ONE HOUR OF COLLECTION FOR PATIENTS ON HEPARIN
PTT CORRECTION STUDY	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE). TUBE MUST BE FULL
PTT INSENSITIVE	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
PYRUVATE DEHYDROGENASE COMPLEX	CONSULT BIOCHEMIST BEFORE PROCUREMENT
PYRUVATE KINASE	AT LEAST 1.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) COLLECT BEFORE 1100 HRS MONDAY TO WEDNESDAY
Q FEVER AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROVIDE CLINICAL HISTORY OF PATIENT
QUANTIFERON-TB	CALL PUBLIC HEALTH LAB FOR TUBES KEEP AT ROOM TEMPERATURE AND DELIVER TO LAB IMMEDIATELY
QUINIDINE	COLLECT 3.5 ML OF BLOOD IN SST TUBE. COLLECT SPECIMEN PRIOR TO LAST DOSE
RABIES SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
RAPAMYCIN (SIROLIMUS)	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
RBC CAST	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
RBC FOLIC ACID ORDER CBC	COLLECT 4.0 ML OF BLOOD EDTA TUBE (LAVENDER) ALSO
RED CELL CHOLINESTERASE	COLLECT 2 X 4.0 ML OF BLOOD EDTA TUBE (LAVENDER)
RENIN	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) CENTRIFUGE AND FREEZE WITHIN 4 HOURS
REPTILASE TIME	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
RESPIRATORY VIRUS	NASOPHARYNGEAL SWAB OR ASPIRATE IN UTM TRANSPORT MEDIUM: INCLUDES INFLUENZA, PARAINFLUENZA, RSV AND ADENOVIRUS
RHEUMATOID FACTOR	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
RICKETTSIAL SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
RIVOTRIL (CLONAZEPAM)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
RSVR	COLLECT NASOPHARYNGEAL SWAB OR ASPIRATE IN UTM TRANSPORT MEDIUM
RUBELLA IgG	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
RUBELLA IgM	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
SACCHAROMYCES CER. AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
SALICYLATE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE TIME, DOSAGE AND TIME COLLECTED
SALIVARY CORTISOL	CALL 777-2373 FOR PROPER COLLECTION TUBE IN ADVANCE
SCALPH PH	COLLECT ON ONE FULL CAPILLARY TUBE THIS IS A PART OF SCALP PH KIT
SCHISTOSOMA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROVIDE CLINICAL HISTORY OF PATIENT
SEDIMENTATION RATE	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
SELENIUM	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
SEMEN ANALYSIS	COLLECT SEMEN SAMPLE IN STERILE URINE CONTAINER. ABSTAIN FROM EJACULATION FOR 3 DAYS KEEP AT 37°C  DELIVER ASAP. APPOINTMENT NECESSARY. DELIVER TO FERTILITY LAB 35 MAJOR'S PATH ST. JOHN'S OR YOUR  LOCAL HOSPITAL IN RURAL NL.
SEROTONIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE  SPECIAL DIET REQUIRED CALL 777-2273
SERUM FERRITIN	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
SERUM FREE CARNITINE	COLLECT 3 LITHIUM HEPARIN CONES (LIGHT GREEN) OR 3.5 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
SERUM FRUCTOSE	COLLECT SEMEN SAMPLE IN STERILE URINE CONTAINER. ABSTAIN FROM EJACULATION FOR 3 DAYS KEEP AT 37°C  DELIVER ASAP APPOINTMENT NECESSARY. DELIVER TO FERTILITY LAB 35 MAJOR'S PATH ST. JOHN'S OR YOUR  LOCAL HOSPITAL IN RURAL NL
SERUM PROTEIN ELECTROPHORESIS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
SG (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR 4.0 ML OF BLOOD IN RED STOPPER TUBE
SICKLE CELL SCREEN	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) PLEASE DO NOT COLLECT ON FRIDAY'S
SILICA CLOTTING TIME	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
SILVER	COLLECT 4 ML OF BLOOD IN SODIUM HEPARIN TUBE (DARK GREEN)
SMOOTH MUSCLE AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
SODIUM	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
SODIUM (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIM TUBE (LIGHT GREEN)
SOL TRANSFERRIN RECEPTOR	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
SOMATOSTATIN	TEST BY APPOINTMENT ONLY. PLEASE CALL 777-7750 OR 777-6418
STONE ANALYSIS	COLLECT STONES IN A CONTAINER SEND TO RENAL LAB
STOOL FOR PH	COLLECT FRESH RANDOM STOOL IN STERILE CONTAINER. STOOL "MUST" BE WATERY CONSISTENCY
STOOL REDUCING SUBSTANCES	COLLECT FRESH RANDOM STOOL IN STERILE CONTAINER FREEZE SPECIMEN IF DELIVERY DELAYED
STRIATED MUSCLE AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
STRONGYLOIDES SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROVIDE CLINICAL HISTORY OF PATIENT
SUCROSE LYSIS	ARRANGE WITH HAEMATOLOGY LAB FOR TESTING INSTRUCTIONS AT 777-6386
SWEAT TEST	ARRANGE WITH MASS SPECTROMETRY LAB FOR TESTING AT 777-4513
SYPHILIS TP AB	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
T1 (MAPLE TREE)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE



**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
T3 (BIRCH)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
TABC	COLLECT 2X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
TAENIA SEROLOGY	COLLECT 3.5 ML OF BLOOD SST TUBE (GOLD) PROVIDE CLINICAL HISTORY OF PATIENT
TB CULTURE	SPUTUM, URINE, BRONCHAIL WASH/LAVAGE, TISSUE, FLUID ETC. IN STERILE CONTAINER. INCLUDES: CULTURE AND AURAMINE SMEAR
TESTOSTERONE TOTAL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
THEOPHYLLINE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT AT TROUGH OR PEAK NOTE TIME, DOSAGE AND TIME COLLECTED
THIOCYANATE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: INCLUDE TIME OF LAST DOSAGE AND COLLECTION TIME
THIOPURINE METHYLTRANSFERASE	COLLECT 3.5 ML OF BLOOD IN EDTA TUBE (LAVENDER) ONLY COLLECT MONDAY, TUESDAY AND BEFORE 10:00 AM ON WEDNESDAY
THROMBIN TIME	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
THYROGLOBULIN+THYROGLOBULIN AB	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
THYROID STIMUL.HORMONE	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) SST TUBE (GOLD) FOR NEONATES TAKE 2 FULL CONES

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
THYROXINE BINDING GLOBULI	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TISSUE ARSENIC	COLLECT 1.0 g OF HAIR AND OR NAIL CLIPPINGS
TISSUE TRANSGLUTAMINASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT ON ICE AND DELIVER TO LAB ASAP REPLACING ANTI-ENDOMYSIAL ANTIBODIES
TISSUE TYPING BMT	COLLECT 2X 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER)
TOBRAMYCIN PEAK	COLLECT 3.5 ML OF BLOOD IN SST TUBE. (GOLD) COLLECT 30 MINUTES AFTER INFUSION NOTE: RECORD TIME OF COLLECTION
TOBRAMYCIN TROUGH	COLLECT 3.5 ML OF BLOOD IN SST TUBE. (GOLD) COLLECT PRIOR TO DOSE. NOTE: RECORD TIME OF COLLECTION
TOLBUTAMIDE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN)
TOPIRAMATE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TOTAL BILI FLUID	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR RED STOPPER TUBE
TOTAL BILIRUBIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TOTAL PROTEIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TOTAL URINE ESTROGEN	COLLECT URINE IN A 24 HR URINE CONTAINER STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
TOXOCARA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROVIDE CLINICAL HISTORY OF PATIENT
TOXOPLASMA PCR	CALL MICROBIOLOGY LAB FOR INSTRUCTIONS AT 777-6936 SAMPLE MUST BE REFERRED TO MONTREAL <24 HOURS
TOXOPLASMA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TOXOPLASMA SEROLOGY IgG	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TOXOPLASMA SEROLOGY IgM	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TRANSFERRIN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TRAZODONE (DESYREL)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TRICHINELLA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROVIDE CLINICAL HISTORY OF PATIENT
TRIGLYCERIDE (FLUID)	COLLECT 3.0 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR RED STOPPER TUBE
TRIGLYCERIDES	COLLECT 3.5 ML OF BLOOD IN SST TUBE 12 HR. FAST REQUIRED
TRIIODOTHYRONINE: TOTAL	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TRIMIPRAMINE	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM (DARK GREEN)
TROPONIN	COLLECT 3.0 ML BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN)
TRYPSINOGEN	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
TRYPTASE	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TSH RECEPTOR AUTOANTIBODIES	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TUMOR LYSIS SYNDROME	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
TYROSINE FILTER PAPER	COLLECT BLOOD ON NEWBORN SCREENING FILTER PAPER
U.COCAINE CONFIRMATION	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
UR.BILIRUBIN	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
UR OPIATE CONFIRMATION	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
UR.AMPHETAMINE	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
UR.BARBITURATE	COLLECT 20 ML OF RANDOM URINE IN STERILE URINE CONTAINER
UR.BENZODIAZEPINES	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
UR.CANNABINOIDS	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
UR.COCAINE	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
UREA	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
UREA (FLUID)	COLLECT 3.5 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR RED STOPPER TUBE

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
UREAPLASMA	URETHRAL/UROGENITAL/VAGINAL/RECTAL SWAB OR URINE IN UTM. NEONATAL BRONCH LAVAGE, THROAT/NP SWAB  OR NP ASPIRATE  RELEVANT HISTORY REQUIRED
URIC ACID (FLUID)	COLLECT 3.5 ML OF FLUID IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR RED STOPPER TUBE
URINALYSIS & MICROSCOPIC	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINARY COPPER	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER  ADJUST 100 ML ALIQUOT PH TO 2-3 WITH 6N HCL AFTER COLLECTION IF SAMPLE IS BEING REFERRED IN FROM OUTSIDE ST. JOHN'S.
URINARY FREE CORTISOL	URINE TO BE COLLECTED IN A 24 HR. PRE-ACIDIFIED CONTAINER (50% ACETIC ACID)  NOTE: PLEASE STATE VOLUME AND TIME OF COLLECTION ON CONTAINER
URINARY ROHPNOL	COLLECT 5 ML OF RANDOM URINE IN NON-STERILE URINE CONTAINER
URINE 5-HIAA	COLLECT URINE IN A PRE-ACIDIFIED (50% ACETIC ACID) 24 HR. URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
URINE ALDOSTERONE	URINE TO BE COLLECTED IN A 24 HR. PRE-ACIDIFIED CONTAINER (50% ACETIC ACID)  NOTE: PLEASE STATE VOLUME AND TIME OF COLLECTION ON CONTAINER
URINE AMINOACIDS QUANTITATIVE	COLLECT 20 ML OF RANDOM URINE IN STERILE CONTAINER

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
URINE AMYLASE	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, TIME AND DATE COLLECTED ON CONTAINER. REFRIGERATE
URINE BARBITURATE SCREEN	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE CADMIUM	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
URINE CALCIUM	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
URINE CARNITINE	COLLECT 10 ML OF URINE IN STERILE URINE CONTAINER
URINE CHLORIDE	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
URINE CITRATE	URINE TO BE COLLECTED IN A 24 HR. PRE- ACIDIFIED CONTAINER (50% ACETIC ACID)  NOTE: PLEASE STATE VOLUME AND TIME OF COLLECTION ON CONTAINER
URINE COPROPORPHYRIN	COLLECT URINE IN A 24 HR SODIUM CARBONATE URINE CONTAINER  STATE VOLUME, DATE AND TIME ON CONTAINER  PROTECT FROM LIGHT
URINE URINE CREAT.mmol/d	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
URINE CREAT.mmol/L	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE CRYSTALS	COLLECT 20 ML OF RANDOM URINE IN STERILE CONTAINER

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
URINE DELTA ALA	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER  PROTECT FROM LIGHT
URINE DESIPRAMINE	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE ETHANOL	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE GAMA HYDROXY BUTYRATE	COLLECT 20 ML OF RANDOM URINE IN STERILE URINE CONTAINER  FREEZE ASAP AND DELIVER TO LAB
URINE GLUCOSE	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE GLUCOSE 24H	COLLECT IN 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER. REFRIGERATE & DELIVER TO LAB
URINE GLUCOSE RANDOM	COLLECT 20 ML OF URINE IN A RANDOM URINE CONTAINER  REFRIGERATE & DELIVER TO LAB
URINE HOMOCYSTEINE	COLLECT URINE IN 24 HR URINE CONTAINER
URINE HYDROXYPROLINE	COLLECT URINE IN 24 HR URINE CONTAINER
URINE INORG.PHOS.	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
URINE IODINE	COLLECT URINE IN A 24 HR URINE CONTAINER, RANDOM URINE ALSO ACCEPTABLE  STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
URINE KETONES	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER.

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
URINE LEAD	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
URINE MAGNESIUM	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
URINE MICROALBUMIN	COLLECT URINE IN 24 HR. CONTAINER OR RANDOM URINE IN STERILE URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
URINE MYOGLOBIN	COLLECT 6.0 ML OF RANDOM URINE IN STERILE URINE CONTAINER
URINE NITRITE	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE OCCULT BLOOD	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE OLIGOSACCHARIDES	COLLECT 20 ML OF RANDOM URINE IN NON-STERILE URINE CONTAINER
URINE ORGANIC ACIDS	COLLECT 5 ML OF RAMDOM URINE IN NON-STERILE URINE CONTAINER  FREEZE IMMEDIATELY
URINE OSMOLALITY	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER  PLEASE DO NOT ADD ACID
URINE OXALATE	URINE TO BE COLLECTED IN A 24 HR. PRE-ACIDIFIED CONTANIER (50% ACETIC ACID)  NOTE: PLEASE STATE VOLUME AND TIME OF COLLECTION ON CONTAINER



**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
URINE PBG QUANTITATION U	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER  PROTECT FROM LIGHT
URINE PORPHYRIN SCREEN	COLLECT URINE IN 24 HR SODIUM CARBONATE URINE CONTAINER  PROTECT FROM LIGHT, REFRIGERATE
URINE POTASSIUM	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
URINE PROTEIN	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
URINE PROTEIN CONCENTRATION	COLLECT URINE IN A 24 HR URINE CONTAINER  STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
URINE PROTEIN g/L	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE REDUCING SUBSTANCES	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE SALICYLATE	COLLECT 20 ML OF RANDOM URINE IN STERILE URINE
URINE SEROTONIN	COLLECT IN 24 HR URINE CONTAINER WITH HCL PRESERVATIVE  STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
URINE SODIUM	COLLECT URINE IN A 24 HR URINE CONTAINER
STURINE SULFITE	CALL DR. ED RANDELL AT 777- 6375 FOR INSTRUCTIONS
URINE TRICYCLIC SCREEN	COLLECT 20 ML OF FRESH URINE IN STERILE URINE CONTAINER

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
URINE UREA	COLLECT URINE IN A 24 HR URINE CONTAINER STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
URINE URIC ACID	COLLECT URINE IN A 24 HR URINE CONTAINER STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
URINE UROBILINOGEN	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE UROPORPHYRIN	COLLECT URINE IN A 24 HR SODIUM CARBONATE URINE CONTAINER PROTECT FROM LIGHT. REFRIGERATE
URINE WBC	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
URINE ZINC	COLLECT URINE IN A 24 HR URINE CONTAINER STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER ADJUST 100 ML ALIQUOT PH TO 2-3 WITH 6N HCL AFTER COLLECTION IF SAMPLE IS BEING REFERRED IN
UR.LT.CHAIN ELECT.	COLLECT URINE IN 24 HR. CONTAINER OR RANDOM URINE IN STERILE CONTAINER STATE VOLUME, DATE AND TIME OF COLLECTION ON CONTAINER
UR.MMA	COLLECT 10 ML OF URINE NON-STERILE URINE CONTAINER
UR.OPIATES	COLLECT 20 ML OF RANDOM URINE IN STERILE URINE CONTAINER
UR.PCP	COLLECT 20 ML OF RANDOM URINE IN STERILE URINE CONTAINER
UR.RITALIN SCREEN	COLLECT 20 ML OF RANDOM URINE IN STERILE URINE CONTAINER

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
VALPROIC ACID	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT PRIOR TO NEXT DOSE NOTE: RECORD TIME OF COLLECTION
VANCOMYCIN PEAK	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT 60 MINUTES AFTER COMPLETING INFUSION NOTE: RECORD TIME OF COLLECTION
VANCOMYCIN TROUGH	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) COLLECT PRIOR TO NEXT DOSE NOTE: TIME RECORD OF COLLECTION
VARICELLA ZOSTER	COLLECT VESICULE FLUID OR LESION IN UTM TRANSPORT MEDIUM
VARICELLA ZOSTER IgG SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
VARICELLA ZOSTER IgM SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
VARICELLA ZOSTER PCR	VESICULAR FLUID OR SCRAPING FROM BASE OF LESION OR BIOPSY MATERIAL IN UTM TRANSPORT MEDIUM; CSF IN STERILE CONTAINER
VASOACTIVE INTEST. POLYPEPTID E	COLLECT 2 X 4.0 ML IN EDTA TUBE (LAVENDER). DELIVER TO LAB ON ICE
VENOUS BLOOD GASES	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN (TUBE) DARK GREEN. DELIVER ON ICE TO LAB IMMEDIATELY.
VENOUS CORD GAS	COLLECT 0.5 ML OF BLOOD IN HEPARINIZED SYRINGE DELIVER ON ICE TO LAB IMMEDIATELY

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
VISCOSITY	COLLECT 10 ML OF BLOOD IN PLAIN RED TUBE
VITAMIN A	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROTECT FROM LIGHT
VITAMIN B1	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) PROTECT FROM LIGHT
VITAMIN B12	COLLECT 4.0 ML OF BLOOD IN LITHIUM HEPARIN TUBE (LIGHT GREEN) OR SST TUBE (GOLD)
VITAMIN B2	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN) PROTECT FROM LIGHT
VITAMIN B6	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) PROTECT FROM LIGHT
VITAMIN C	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIUM TUBE (DARK GREEN) PROTECT FROM LIGHT
VITAMIN E	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) PROTECT FROM LIGHT
VITAMIN K	COLLECT 3.0 ML OF SERUM IN SST TUBE (GOLD) PROTECT FROM LIGHT
VITREOUS ETHANOL	COLLECT 3.5 ML OF FLUID IN SODIUM FLORIDE TUBE (GREY) TO BE USED FOR FORENSIC SPECIMENS ONLY
VMA	URINE TO BE COLLECTED IN A 24 HR. PRE-ACIDIFIED CONTAINER (50% ACETIC ACID)  STATE DATE, VOLUME AND TIME OF COLLECTION ON CONTAINER SEND TO LAB IMMEDIATELY

**EASTERN HEALTH TEST DICTIONARY**

<b>TEST NAME</b>	<b>INSTRUCTIONS</b>
VOLUME in mLs	COLLECT URINE IN A 24 HR URINE CONTAINER STATE VOLUME, DATE AND TIME COLLECTED ON CONTAINER
VON WILLEBRANDS GROUP	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
VWF MUTATION	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
W8 (DANDELION)	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD) NOTE: UP TO FIVE ALLERGENS CAN BE DONE ON ONE TUBE
WAXY CAST	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
WBC CAST	COLLECT 20 ML OF URINE IN STERILE URINE CONTAINER
WEST NILE VIRUS	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
XA AB	COLLECT 4.5 ML OF BLOOD IN NA CITRATE TUBE (BLUE) TUBE MUST BE FULL
XYLOCAINE	COLLECT 3.0 ML OF SERUM IN SST (GOLD)
YERSINIA SEROLOGY	COLLECT 3.5 ML OF BLOOD IN SST TUBE (GOLD)
ZINC	COLLECT 4.0 ML OF BLOOD IN SODIUM HEPARIN (DARK GREEN)
ZINC PROTOPHYRIN	COLLECT 4.0 ML OF BLOOD IN EDTA TUBE (LAVENDER) DELIVER TO LAB ASAP